



Deluxe Connected Blood Pressure Monitor UA-651BLE



Introducing the Deluxe Blood Pressure monitor UA-651BLE, one of the most advanced monitors available designed for ease of use and accuracy with a simple, compact and modern design. The UA-651BLE offers advanced features including Pressure Rating Indicator™, irregular heartbeat detection, low battery indicator, measurement time and date stamp and up to 30 readings in internal memory storage.

Unlock the Potential of the UA-651BLE with the A&D WellnessConnected™ Platform

You can now achieve greater insight into your health and wellness like never before. As part of A&D Medical's WellnessConnected™ portfolio, the UA-651BLE easily, reliably and securely connects to your iPhone or Android supported mobile device with the A&D Connect app via Bluetooth® Smart™ technology. A&D Connect is a simple, easy to use mobile application that allows you to view your wellness metrics in a whole new way, complete with trends and graphs for all of your A&D WellnessConnected™ devices, including blood pressure, weight, activity and sleep. With an A&D WellnessConnected™ account, you have unlimited storage in the A&D Cloud.

Track your wellness. Share with others.

A&D Connect allows users to engage physicians and caregivers in their health and wellness monitoring. Set your goals, see your trends and live well. The A&D Connect sharing function allows you to share your health and wellness data with family, friends and caregivers via email, social media platforms and more.



A&D Connect is compatible with Bluetooth® Smart Ready™ mobile devices from Apple and many Android device manufacturers.

APPLE:




iOS version 6 or newer

iPhone 4s or newer, iPad 3rd & 4th generation or newer, iPad Mini

ANDROID:

OS version 4.4 or newer

Requires Bluetooth Smart Ready™ device (Bluetooth 4.0).

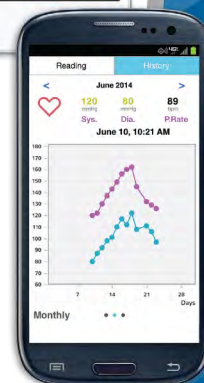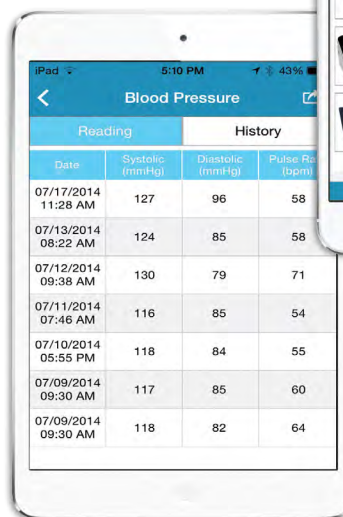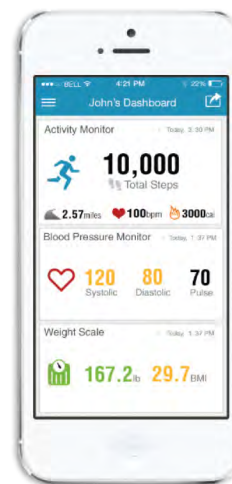

For a list of supported devices, visit www.wellnessconnected.com/mobile



- Clinically Validated — using latest protocol ESH 2010
- SlimFit™ Cuff System w/ Gentle inflation — Flexible inflation gives less chance of re-inflation
- Advanced Measurement and Report
 - 30 reading internal memory storage / Unlimited reading storage with A&D Cloud account
 - Stores memory with time and date stamp
 - Multiple users via A&D Connect app
 - A&D Connect app allows for easy sharing of readings via mobile social applications, email and more.
- Risk Factor Indicators
 - Gives a snap shot of your BP classification in accordance with US JNC7 guidelines.
 - Irregular Heartbeat Feature – Alerts the user to the presence of an irregular heartbeat risk

Part Number	UPC Number	Cuff Size	*Included. Note: Model UA-651BLE is not designed to use a small cuff.
UA-651BLE	093764603327		
Medium SlimFit™ - UA-290*	093764602986	9.0" – 14.6" (23cm – 37cm)	
Large SlimFit™ - UA-291	093764602993	12.2" – 17.7" (31cm – 45cm)	

Included Accessories:

- Four AA batteries
- Bilingual Quick Start Card
- Bilingual Instruction Manual (English and Spanish)
- Carrying case

TECHNICAL SPECIFICATIONS	
Measurement method	Oscillometric measurement
Measurement range	Pressure: 20 - 280 mmHg Systolic pressure: 60 - 279 mmHg Diastolic pressure: 40 - 200 mmHg Pulse: 40 - 180 beats /min
Measurement accuracy	Pressure: ±3 mmHg Pulse: ±5 %
Power supply	4 x 1.5V batteries (R6P, LR6 or AA) or AC adapter (TB:233 - sold separately)
Number of measurements	Approx. 700 times LR6 (alkaline batteries) Approx. 200 times R6P (manganese batteries) With pressure value 180 mmHg, room temperature 23 °C.
Classification	Internally powered ME equipment (Supplied by batteries) / Class II (Supplied by adapter) Continuous operation mode
Clinical test	According to ANSI / AAMI SP-10 1992
Wireless communication	VZ (MURATA Manufacturing Co.Ltd.) Bluetooth Ver. 4.0LE BLP

PACKAGING AND SHIPPING SPECIFICATIONS			
Individual Unit (Box)	W x D x H (inches)	WT (lb.)	Individual Item Cube (cubic feet)
UA-651 BLE	5.12 x 4 x 7.64	1.52	.088
Unit Case	W x D x H (inches)	WT (lb.)	Case Cube (cubic feet)
10 units	21.65 x 11 x 8.46	15.87	1.18